
















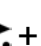



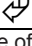


Function	Description	Internal Keyboard (S1...S6 = memory key)	PC-Keyboard	Remarks
Entry	Record text in Display	Standard keys	Standard keys	
Backspace	Delete the character to the left of the cursor		Backspace	
Cursor left	Move the cursor left (back)	 + , (←)	←	For repeat, keep Shift pressed (internal keyboard)
Cursor right	Move cursor right (forward)	 + . (→)	→	For repeat, keep Shift pressed (internal keyboard)
Home	To the start of text	 + , (←)	Home	
End	To the end of text	 + . (→)	End	
Clear the display	All text in the display is cleared	 + space	Alt + Space or F8	
Shift lock	Set capital letters	2 x 	Press shift twice	LED flashes, then remains lit on Shift Lock
Speak	Speak display contents		F10	
Store message	Store message for speaking immediately	 +  then S1 - S6	Alt + Alt then F2 - F7	Stored in current field
Store sentence	Store sentence for processing later	 +  then number or letter	Alt + Alt then number or letter	Stored in current field
Recall message	Directly speak message saved from level 1	S1 – S6	F2 – F7	Spoken immediately
Load sentence	Load saved sentence to display at cursor position	 + number or letter	Alt + number or letter	Can be spoken by pressing the Speak key
Change field level 1	A new communication field from level 1 is to be selected	 +  (= field menu) then S1 - S6	Alt + F10 (= field menu) then F2 - F7	6 fields appear in the display, each next to the memory keys. By pressing the key the respective field is selected.
Change field level 2	A new communication field from level 2 is to be selected	 +  +  (= field menu) then S1 - S6	Shift + Alt + F10 (= field menu) then F2 - F7	Another 6 fields appear in the display, each next to the memory keys. By pressing the key the respective field is selected.
Call setups menu		 + 	Alt + Enter	5 menu points and 'back' appears in the display
Signal sound			F9	
End settings	Return from any menu point		Enter	
Switch off unit		Keep one of the 6 direct memories pressed down		Until the text 'Release key' appears